

ADscope™ Cardiology



Adscope™ 606 Ultra-Light Cardiology Stethoscope

Cardiology scope performance at a fraction of the weight.

The 606 series features

- 40% lighter than our best selling 601 series Cardiology Scope
- Combination diaphragm/bell chestpiece precisely CNC machined to exacting tolerances from ultra-lightweight aluminum alloy
- Extra-deep bell for unsurpassed low-frequency response
- Ultra-sensitive diaphragm for greater amplification and crisper high-frequency transmission
- Cardiology headset with large bore aluminum binaurals are fixed at 15° angle to maximize comfort and acoustic seal
- Bi-lumen (dual channel) extra-thick 19" PVC tubing optimizes acoustic transmission
- Adsoft™ threaded PVC eartips for the ultimate in wearing comfort and acoustic seal
- Includes accessory kit with spare diaphragm, extra Adsoft™ and standard eartips
- Scope ID Tag included
- Weighs 128 grams (4.5 oz.)
- Overall length 71 cm (28")
- Lifetime Warranty

Ordering Information


| Item | Description |
|-------|--|
| 606** | Adscope™ 606, Cardio-Lite Scope (Specify Color) |

Color Choices



Adscope™ 600 Platinum Cardiology Stethoscope

The sculpted chestpiece and AFD Technology diaphragm provide wide frequency performance in a compact profile.

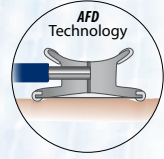

The 600 series features

- Adjustable Frequency Design (AFD) Technology combines the acoustic response of a traditional bell and diaphragm into a convenient one-sided chestpiece
- Light pressure for low bell-type frequency, firm pressure for diaphragm type high frequency
- Sculpted chestpiece precisely CNC machined to exacting tolerances from surgical stainless steel for outstanding performance and rugged durability
- Non-chill diaphragm retaining rim for patient comfort
- Cardiology headset with large bore stainless steel binaurals are fixed 15° angle to maximize comfort and acoustic seal
- Bi-lumen (dual channel) extra-thick 19" PVC tubing optimizes acoustic transmission
- Adsoft™ threaded PVC eartips for the ultimate in wearing comfort and acoustic seal
- Accessory kit includes 1 additional pair of Adsoft™ and standard eartips
- Scope ID Tag included
- Weighs 23 gr (7.15 oz)
- Overall length 68.5cm (27")
- Lifetime Warranty

Color Choices for 600** Only




Ordering Information

| Item | Description |
|--------------|--|
| 600** | Adscope™ 600, Cardiology, (Specify Color) |
| 600GP | Adscope™ 600GP, Cardiology, (Black Tubing Only) |
| 600ST | Adscope™ 600, Cardiology, Tactical (Black Tubing Only) |

Adscope™ 601 Convertible Cardiology Stethoscope

Combining AFD Technology with the convertible design makes this the most versatile of our Cardiology Scopes

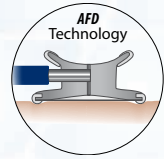

The 601 series features

- Threaded chestpiece drum allows rapid conversion from a pediatric diaphragm to an adult bell
- Adjustable Frequency Design (AFD) Technology on both adult and pediatric diaphragm enhances frequency response
- Chestpiece precisely CNC machined to exacting tolerances from surgical stainless steel for outstanding performance and rugged durability
- Extra deep bell for unsurpassed low frequency response
- Non-chill bell and diaphragm retaining rims for patient comfort
- Cardiology headset with large bore stainless steel binaurals are fixed at 15° angle to maximize comfort and acoustic seal
- Bi-lumen (dual channel) extra-thick 19" PVC tubing optimizes acoustic transmission
- Adsoft™ threaded PVC eartips for the ultimate in wearing comfort and acoustic seal
- Includes accessory kit with bell chestpiece, extra Adsoft™ and standard eartips
- Scope ID Tag included
- Weighs 215 gr (7.6 oz.)
- Overall length 71 cm (28")
- Lifetime Warranty

Color Choices



Ordering Information

| Item | Description |
|--------------|--|
| 601** | Adscope™ 601, Cardiology Scope (Specify Color) |