

The first modular multiplatform diagnostic station that grows with your practice

At blood pressure technology to create the first truly modular diagnostic station that grows with your needs. Start with basic automated BP, then add thermometry and pulse oximetry AS you need it. Choose from wall mount, table mount or mobile configurations. Or convert the unit to the platform needed as your requirements change.

The complete ADVIEW™ diagnostic station includes:

- **An automated BP monitor** that measures systolic and diastolic pressure, pulse rate and MAP* in about 30 seconds using the oscillometric method.
- Our unique manual mode allows measurement with your stethoscope using the auscultatory method.
- Integrated rechargeable battery enhances mobility. Provides up to 200 measurements on a single charge.
- **A temperature module using Filac™** FasTemp™ technology that provides an oral or axillary temperature in about 10 seconds. Built in storage well for the probe and box of 20 disposable probe covers.
- **An SpO₂ module** that features industry leading Nellcor compatible technology for accurate readings in about 10 seconds.
- **Optional EMR** (electronic medical records) with Bluetooth connectivity for wireless integration to your computer system.*

Modular

Unlike EVERY other monitor on the market, the ADVIEW™ is completely modular. You can purchase the ADVIEW™ configured as you want from the factory OR upgrade the unit in the field at any time. No tools needed. With connect-and-go capability, adding optional temperature and pulse oximetry modules is effortless.

Adaptable

With three platforms - wall, tabletop, or mobile - you can choose the modality that's right for you.

* Bluetooth equipped device is EMR ready. Consult with your EMR provider for integration costs and details.



Ordering Information

Adviv Monitor

- * **9000BP** Adviv BP
- * **9000BPPT** Adviv BP/Temp
- * **9000BPS** Adviv BP/SpO₂
- * **9000BPSTO** Adviv BP/SpO₂/Temp

Adviv Monitor with Bluetooth Connectivity

- * **9001BP** Adviv BP with Bluetooth
- * **9001BPPT** Adviv BP/Temp with Bluetooth
- * **9001BPS** Adviv BP/SpO₂ with Bluetooth
- * **9001BPSTO** Adviv BP/SpO₂/Temp with Bluetooth

BP monitor includes: Rechargeable battery; 9v AC adapter with US power cord; 8 ft length of tubing; adult cuff; large adult cuff; Operating manual on CD; Quick Start guide; and ready reference card

Temp module includes: Oral/axillary probe (blue), box of 20 probe covers

SpO₂ module includes: Finger sensor with 10' cable

Complete instrument includes all of the above

Adviv Platforms (Available Separately)

- * **9000D** Desktop caddy
- * **9000W** Wall Mount with Basket
- * **9000M** Mobile Stand

Note:

- *Desk and Wall Mounts shipped inside box with Adviv.
- *Mobile stand shipped separately.

Adviv Modules (Sold Separately)

- * **9000TO** Temperature Module, Oral
- * **9000TR** Temperature Module, Rectal
- * **9000S** SpO₂ Module

* Indicates most popular models



ADC
55 Commerce Drive
Hauppauge, NY 11788

Made in China
Assembled in the U.S.A.

Suntech Medical Ltd.
Oakfield Industrial Estate
Eynsham, Oxfordshire
OX29 4TS UK

tel: 631-273-9600
toll free: 1-800-232-2670
fax: 631-273-9659

www.adctoday.com/adviv
info@adctoday.com

ADVIEWLIT
Printed in U.S.A.



Adviv Cuff Options

- * **850-9000-11AN** Adult cuff, Navy
- * **850-9000-11LAN** Long Adult cuff, Navy
- * **850-9000-9CGR** Child Cuff, Green
- * **850-9000-9LCGR** Long Child Cuff, Green
- * **850-9000-10SARB** Sm. Adult Cuff, Royal
- * **850-9000-10LSARB** Long Sm. Adult Cuff, Royal
- * **850-9000-12XBD** Lg. Adult Cuff, Burgundy
- * **850-9000-12LXBD** Long Lg. Adult Cuff, Burgundy
- * **850-9000-13TBR** Thigh Cuff, Brown
- * **850-9000-13LTBR** Long Thigh Cuff, Brown
- * **9000PCK** Pediatric Cuff Kit (child, sm. adult, adult)
- * **9000ACK** Adult Cuff Kit (adult, lg. adult, thigh)
- * **9000CK** Complete Cuff Kit (child, sm. adult, adult, lg. adult, thigh)

Cuff Ranges

Child	Small Adult	Adult	Large Adult	Thigh
12-19cm	17-25cm	23-33cm	31-40cm	38-50cm

All cuffs equipped with quick lock connector

Adviv Accessories

- * **9000TP** Therm Probe Covers (20/box, 25 boxes/Case)
- * **9000TP-5000** Therm Probe Covers (5000/Case)
- * **9000TR-01** Rectal Temp Kit
- * **9000TRP** Rectal Probe Assembly
- * **9000TO-01** Oral Temp Kit
- * **9000TOP** Oral Probe Assembly
- * **9000SP10** SpO₂ Finger Sensor with 10' cable
- * **9000SPDYS** SpO₂ Pediatric/Multi-Site Sensor (ages 3yrs+)
- * **9000BAT** Rechargeable Battery
- * **9000AC9V** 9v Power Supply
- * **9000PC** Power cord, USA, 120VAC
- * **9000PCEU** Power cord, EURO, 240VAC
- * **9000PCUK** Power cord, UK, 240VAC
- * **9000PCAU** Power cord, AUST, 120VAC
- * **9000QCF-10** Quick Connector, Female, 10/pkg
- * **9000QCM-10** Quick Connector, Male, 10/pkg

Other parts and accessories available

**MAP Available outside U.S.A. only.

For detailed specs visit
www.adctoday.com/adviv
and follow the links

Distributor

INTRODUCING THE NEW

ADVIEW

Blood Pressure
Temperature
Pulse Rate
SpO₂
MAP

Modular, Multiplatform
Affordable Diagnostic Station

Quick, accurate vital signs
at a price that finally
makes sense.



First, Select the CONFIGURATION You Need...

Innovative and easy to use

Purchase the Adview™ configured as you want, and upgrade in the field at any time. No tools required! With our connect-and-go capability, adding optional temperature and pulse oximetry modules is effortless.

Modular Diagnostic Station Base Unit

Blood Pressure:

Non-invasive BP using the oscillometric method. Pulse rate and mean arterial blood pressure (MAP*) also displayed in about 30 seconds. Works with cuff sizes from child to thigh (pediatric patients 3 years and older). Comes equipped with adult and large adult cuffs standard. Other sizes available.

Manual Mode:

For situations in which automatic mode is inappropriate, our unique manual mode allows the operator to measure BP with a stethoscope using the traditional auscultatory method.

*Available outside U.S.A. only.

Optional Modules/Features

Temperature:

A temperature module using Kendall's Filac™ FasTemp™ technology that provides an oral or axillary temperature in about 10 seconds. Scale switchable between Fahrenheit and Centigrade. Built in color coded storage well for the probe and box of 20 disposable probe covers. Easily switched between oral and axillary measurement with the push of a button. Optional color coded rectal module available as well.

Pulse Oximetry:

An SpO₂ module that features industry leading Nellcor compatible technology for accurate readings in about 10 seconds. Complete with Nellcor compatible probe and 10 foot cable. Displays pulse rate signal strength during measurement. Optional pediatric / multi-site sensor available.


EMR Ready Option*:

Upgrade the BP module at time of purchase to include long range (30 meter) Bluetooth connectivity  to your EMR system. Or upgrade later for a one time charge.

*Consult your EMR system provider for integration costs and details.



Adview™ Modular Diagnostic Station
The complete Adview™ monitor incorporates all 3 modules:
Blood Pressure, Temperature, and SpO₂ ...
ALL in one affordable multiparameter device.

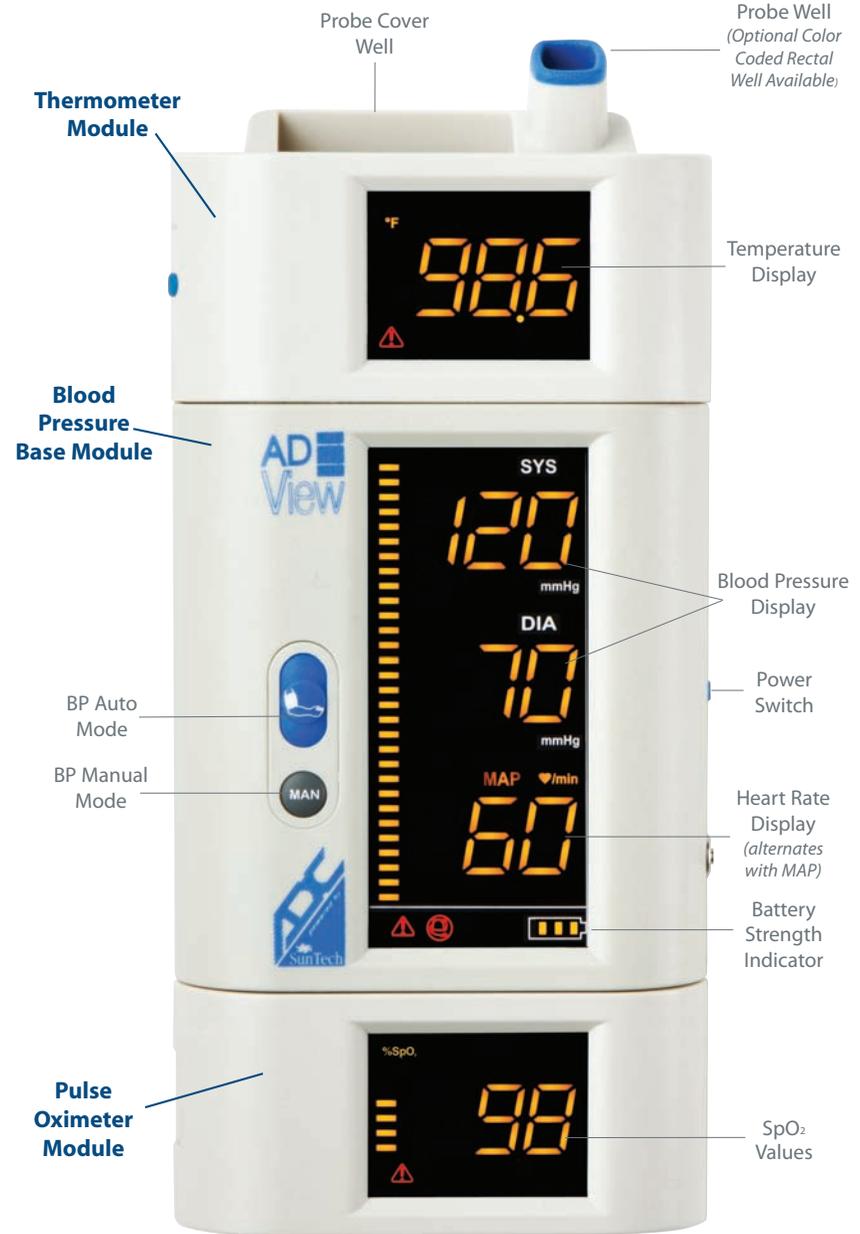

...Then Choose Your PLATFORM

Choose the platform that best suits your needs



Wall Mount

With custom mounting plate that configures quickly and easily to many of the leading brands of clock aneroid. Cuff basket included.



Mobile

With our rugged 5 leg mobile stand you can move the instrument from room to room for optimum workflow.

Table Top

Innovative caddy that holds the unit, power supply and cuff.